

NEW XRT 4000 **Multi Imaging Table**

This trolley is electronically operated and constructed from lightweight material designed to withstand the rigours of hospital life.

Specifically designed for use in Imaging and A&E Departments

The trolley can accommodate conventional cassettes and digital technology, ensuring that new trolleys are not required as imaging equipment is upgraded. Simple detachment of the bucky shelf will allow C arm use.

The patient is able to remain on one trolley during transportation, treatment and X-ray. This saves time, cost and minimises the risk of patient and staff injury as the need to transfer the patient from one trolley to another is avoided.



The XRT4000 has the following features as standard:

- ▶ Maximum weight capacity of 177Kg/390lbs (28 stone)
- ▶ Can be used with conventional, digital and C-arm equipment
- ▶ Trolley is suitable for conventional X-ray cassettes and mobile detectors/plates
- ▶ X-ray facility for the whole length of the table top, maximising the patient's imaging potential
- ▶ The bed raises and lowers via two electric columns driven by a 24V DC battery

The trolley is available in both Adult and Paediatric (XRT4000/P) versions.




XRT 4000 Multi Imaging Table

The full range of Wardray Premise patient transportation equipment can be viewed on our website:

www.wardray-premise.com

- ▶ Supplied with 1 handset and 2 batteries for the electrical operation of the trolley and a charging unit
- ▶ Minimum table height of 725mm, to a maximum of 1025mm
- ▶ Trendelenburg and reverse Trendelenburg

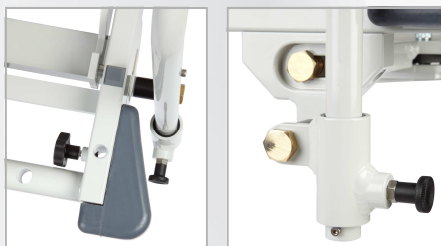

- ▶ Four IV socket positions
- ▶ Drop down cot sides
- ▶ 150mm castors, operated by means of foot bars, with directional lock



- ▶ 50mm mattress with welded seams, to aid infection control procedures
- ▶ The trolley has a 6 position manually operated back rest



Optional extras:

- ▶ Additional handsets for electrical operation of the table
- ▶ Patient's belongings tray (note: restricts C-arm access)
- ▶ Oxygen cylinder holder (note: restricts C-arm access)
- ▶ Foot or head mounted monitor stand
- ▶ Other type of mattress available including memory foam and pressure relief
- ▶ Customised paintwork from our stock RAL colour range
- ▶ Power assisted backrest (note: this will reduce the amount of length that can be X-rayed)

Dimensions:

Length: 2160mm

Width Overall: 700mm

Table Width: 575mm

Table Minimum Height: 725mm

Table Maximum Height: 1025mm

Paediatric Imaging Table XRT4000/P

This trolley has all the features listed for the adult version of the Multi Imaging Trolley with the addition of padded cot side bumpers as standard, to assist with the safe transfer of paediatric patients.

Optional extras:

- ▶ As adult version



Dimensions:

Length: 1500mm

Width Overall: 700mm

Table Width: 575mm

Table Minimum Height: 725mm

Table Maximum Height: 1025mm

GEN.S.14/267A